

▶ Test Menu

▶ Thyroid

- TSH
- Total T3
- Total T4
- FT3
- Ft4

▶ Fertility

- AMH
- Progesterone
- β-HCG
- INH-B*
- PROL
- Testosterone
- FSH
- LH
- E2
- E3*

▶ Tumor Marker

- CA125
- CEA
- tPSA
- fPSA
- AFP
- CA19-9*
- CA15-3*
- PG I
- PG II
- Gastrin-17

▶ Inflammation

- PCT
- CRP
- SAA
- IL-6
- HBP*

▶ Anemia

- VB12
- Folate*
- Ferritin

▶ Cardiac Markers

- cTnI
- Myo
- CK-MB
- NT-ProBNP
- D-Dimer
- cTnT
- H-FABP
- Hs-cTnI
- BNP
- Lp-PLA2*
- St2*
- MPO*
- Hcy*
- GDF-15*

▶ TORCH

- TOXO IgG*
- TOXO IgM*
- RUB IgG*
- RUB IgM*

▶ Diabetes

- Insulin*
- C-Peptide*

▶ Brain Injury

- S100*
- Aβ1-42*
- AD7C-NTP*

▶ COVID-19

- IgG*
- IgM*
- Antigen*
- Neutralizing Antibody (nAb) *

▶ Infection Disease

- TB-IGRA
- HBs-Ag*
- Anti-HCV*
- HIV*
- TP*

▶ Autoimmune

- Anti-TPO*
- Anti-TG*

▶ Bone Metabolism

- 25-OH-VD
- BGP*
- PTH*

* Upcoming

GE BIOMED TECHNOLOGY CO.,LTD

📍 Incheon city, Namdong gu Gansok 3 dong, Gansok L.H.1 danji, 105 dong
130,1 South Korea

☎ +82(0)1067606469

✉ info@ge-biomed.co.kr

🌐 www.ge-biomed.co.kr



Easy Life Plus

Chemiluminescence Immunoassay Analyzer

► EasyLife Plus


Principle	ALP-AMDDP chemiluminescence immunoassay
Sample type	Whole blood/Serum/Plasma
Throughput	30 t/hour
Test mode	Batch mode
Display	10.1-inch touch screen
Calibration frequency	100 days
Consumable	Disposable tips
Data storage	≥100,000
Communication port	USB/RS232/LAN/WIFI, available for LIS
Printer	Built-in thermal printer
Dimensions	575*325*605 mm
Weight	40 kg

► Easy Operation



Step 1

Prepare reagent cartridges
Put reagent cartridges into the sample tray

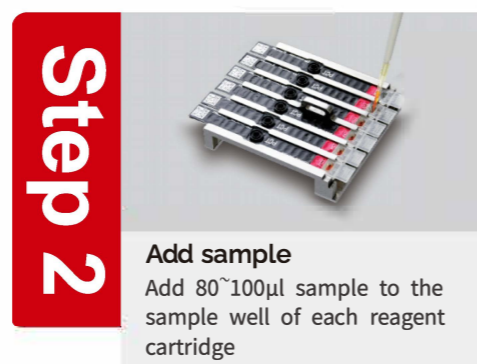


Step 2

Add sample
Add 80~100μl sample to the sample well of each reagent cartridge



Step 3


Insert reagent cartridges
Insert up to 6 reagent cartridges into the analyzer, add one pipette tip for each test

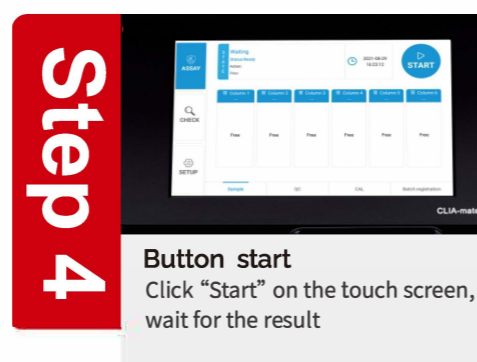

Step 4

Button start
Click "Start" on the touch screen, wait for the result

► All-in-one Reagent Cartridge

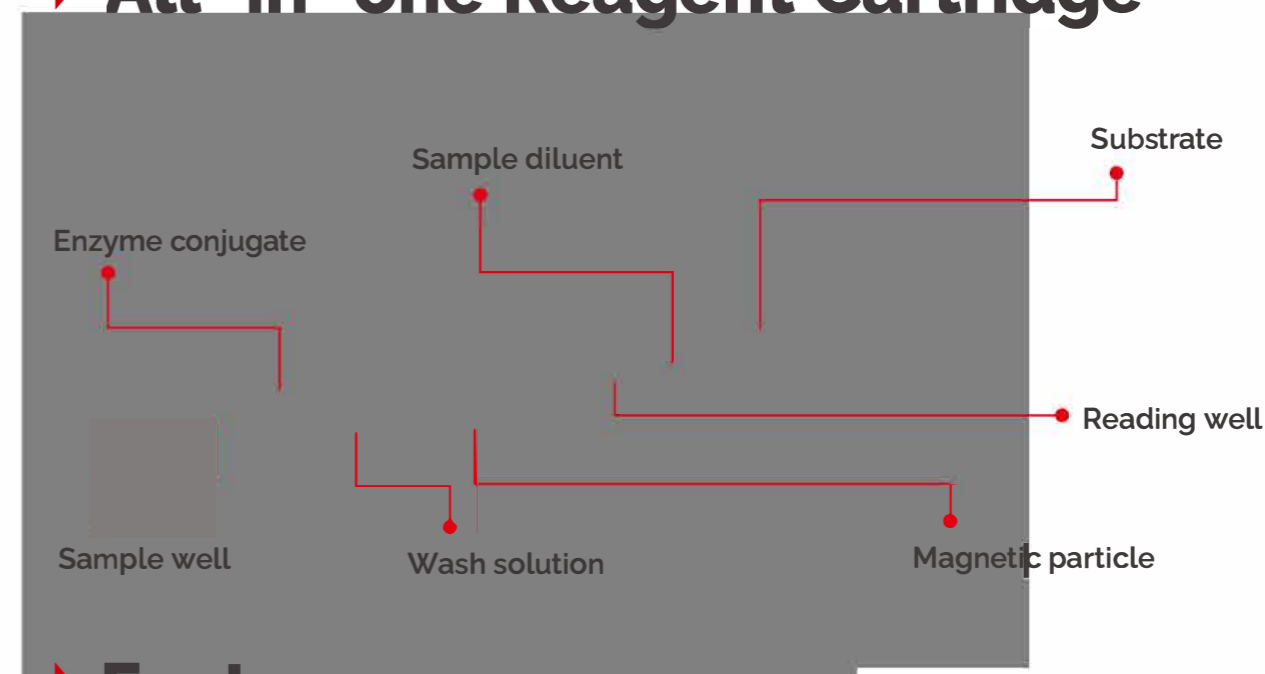

► Feature

- ◎ Real CLIA technology with high sensitivity&accuracy
- ◎ Instrument with compact design at affordable cost
- ◎ Single test reagent cartridge without any reagent waste
- ◎ Specification: 30 test/kit, 60 test/kit
- ◎ Disposable tip for each assay, reduce cross contamination
- ◎ Reagent information loading by external barcode scanner automatically
- ◎ 2 points calibration only, less calibration cost
- ◎ Android operation system, easy upgrading through USB port
- ◎ Easy operation and free maintenance
- ◎ Whole blood for cardiac and inflammation parameters, results within 17 minutes.



► Principle

